

University of Groningen

Suitability of the Cellient (TM) cell block method for diagnosing soft tissue and bone tumors

Song, W.; van Hemel, B. M.; Suurmeijer, A. J. H.

Published in:
Diagnostic Cytopathology

DOI:
[10.1002/dc.23887](https://doi.org/10.1002/dc.23887)

IMPORTANT NOTE: You are advised to consult the publisher's version (publisher's PDF) if you wish to cite from it. Please check the document version below.

Document Version
Publisher's PDF, also known as Version of record

Publication date:
2018

[Link to publication in University of Groningen/UMCG research database](#)

Citation for published version (APA):

Song, W., van Hemel, B. M., & Suurmeijer, A. J. H. (2018). Suitability of the Cellient (TM) cell block method for diagnosing soft tissue and bone tumors. *Diagnostic Cytopathology*, 46(4), 299-305.
<https://doi.org/10.1002/dc.23887>

Copyright

Other than for strictly personal use, it is not permitted to download or to forward/distribute the text or part of it without the consent of the author(s) and/or copyright holder(s), unless the work is under an open content license (like Creative Commons).

The publication may also be distributed here under the terms of Article 25fa of the Dutch Copyright Act, indicated by the "Taverne" license. More information can be found on the University of Groningen website: <https://www.rug.nl/library/open-access/self-archiving-pure/taverne-amendment>.

Take-down policy

If you believe that this document breaches copyright please contact us providing details, and we will remove access to the work immediately and investigate your claim.

Downloaded from the University of Groningen/UMCG research database (Pure): <http://www.rug.nl/research/portal>. For technical reasons the number of authors shown on this cover page is limited to 10 maximum.

ORIGINAL ARTICLE

Suitability of the Cellient™ cell block method for diagnosing soft tissue and bone tumors

W. Song MD | B. M. van Hemel MD | A. J. H. Suurmeijer MD, PhD 

Department of Pathology and Medical Biology, University Medical Center Groningen, University of Groningen, P.O. Box 30.001, RB Groningen, 9700, The Netherlands

Correspondence

A. J. H. Suurmeijer, Department of Pathology and Medical Biology, University Medical Center Groningen, University of Groningen, P.O. Box 30.001, 9700RB Groningen, The Netherlands.
E-mail: a.j.h.suurmeijer@umcg.nl

Funding information

China Scholarship Council (CSC) program, Grant/Award Number: 201606940023

BACKGROUND: The diagnosis of tumors of soft tissue and bone (STB) heavily relies on histological biopsies, whereas cytology is not widely used. Cellient™ cell blocks often contain small tissue fragments. In addition to Hematoxylin and Eosin (H&E) interpretation of histological features, immunohistochemistry (IHC) can be applied after optimization of protocols. The objective of this retrospective study was to see whether this cytological technique allowed us to make a precise diagnosis of STB tumors.

METHODS: Our study cohort consisted of 20 consecutive STB tumors, 9 fine-needle aspiration (FNAC) samples, and 11 endoscopic ultrasonography (EUS) FNACs and included 8 primary tumors and 12 recurrences or metastases of known STB tumors.

RESULTS: In all 20 cases, H&E stained sections revealed that diagnostically relevant histological and cytological features could be examined properly. In the group of 8 primary tumors, IHC performed on Cellient™ material provided clinically important information in all cases. For instance, gastrointestinal stromal tumor (GIST) was positive for CD117 and DOG-1 and a PEComa showed positive IHC for actin, desmin, and HMB-45. In the group of 12 secondary tumors, SATB2 was visualized in metastatic osteosarcoma, whereas expression of S-100 was present in 2 secondary chondrosarcomas. Metastatic chordoma could be confirmed by brachyury expression. Two metastatic alveolar rhabdomyosarcomas were myf4 positive, a metastasis of a gynecologic leiomyosarcoma was positive for actin and estrogen receptor (ER) and a recurrent dermatofibrosarcoma protuberans expressed CD34.

CONCLUSION: In the proper clinical context, including clinical presentation with imaging studies, the Cellient™ cell block technique has great potential for the diagnosis of STB tumors.

KEYWORDS

cell block, Cellient, cytopathology, immunochemistry, soft tissue and bone tumors

1 | INTRODUCTION

Soft tissue and bone (STB) tumors are rare malignancies, which comprise approximately 2% of all neoplasms. Because of this low incidence, patients with STB tumors are usually referred to expert sarcoma centers, where multidisciplinary teams, according to well-established protocols and latest developments, can perform diagnostic procedures and treatment. The clinical diagnosis of primary STB tumors relies on imaging studies and an adequate biopsy. Imaging, in particular with

computed tomography (CT) and magnetic resonance imaging (MRI), provides insight into the location, size, margins, and tissue composition and heterogeneity of STB tumors. Although cytology has been applied in just a few sarcoma centers,^{1–7} the primary diagnosis of STB tumors is usually made on histological (needle or open) biopsies, since these malignancies are morphologically heterogeneous and several histological types have overlapping microscopic features. Moreover, for a conclusive diagnosis of STB tumors, additional immunohistochemistry (IHC) and molecular pathology (fluorescent in situ hybridization (FISH),

This is an open access article under the terms of the Creative Commons Attribution-NonCommercial-NoDerivs License, which permits use and distribution in any medium, provided the original work is properly cited, the use is non-commercial and no modifications or adaptations are made.

© 2018 The Authors. Diagnostic Cytopathology Published by Wiley Periodicals, Inc.

TABLE 1 The 20 STB tumors included in this study

Primary tumor	Localization	Material	Diagnostic antibodies	Diagnosis
1	Stomach	EUS	CD117, DOG1	GIST
2	Stomach	EUS	CD117, DOG1	GIST
3	Peripancreatic	EUS	CD117, DOG1	GIST
4	Stomach	EUS	DOG1	GIST
5	Rectum	EUS	CD117, DOG1	GIST
6	Pararectal	EUS	Beta-catenin	Desmoid fibromatosis
7	Rectum	EUS	Actin, desmin, HMB-45	PEComa
8	Retroperitoneum	EUS	Desmin, caldesmon, SMA	Leiomyosarcoma
Secondary tumor				Metastasis or recurrence of
1	Mediastinum	EUS	Desmin, ER	Uterine leiomyosarcoma
2	Pancreas	EUS	SATB2	Osteosarcoma of bone
3	Mediastinum	EUS	S-100	Chondrosarcoma of bone
4	Inguinal node	FNAC	Brachyury	Chordoma of sacrum bone
5	Inguinal node	FNAC	...	Pleomorphic sarcoma NOS
6	Orbit	FNAC	myf4	Alveolar rhabdomyosarcoma
7	Cheek	FNAC	myf4	Alveolar rhabdomyosarcoma
8	Neck	FNAC	...	Radiation-induced MPNST
9	Inguinal node	FNAC	...	Pleomorphic sarcoma NOS
10	Paranasal	FNAC	CD-34	Dermatofibrosarcoma (DFSP)
11	Hip	FNAC	S-100	Chondrosarcoma of bone
12	Supraclavicular	FNAC	...	Pleomorphic radiation sarcoma

Abbreviations: DFSP, dermatofibrosarcoma protuberans; GIST, gastrointestinal stromal tumor; MPNST, malignant peripheral nerve sheath tumor; NOS, not otherwise specified.

polymerase chain reaction (PCR), and next generation sequencing (NGS)) often have to be administered, requiring special expertise.

Typing and grading of STB tumors is mandatory for treatment decisions. In this clinical context, in our and most other sarcoma teams, cytology is only applied in selected cases. Firstly, in cases with an established diagnosis of the primary STB tumor, cytology can be effectively used to diagnose recurrent or metastatic sarcoma. Secondly, for deep-located STB tumors, cell material can be collected by fine-needle aspiration (FNAC) during endoscopic ultrasonography (EUS-FNA) or by endobronchial ultrasound-guided transbronchial needle aspiration (EBUS). In our sarcoma team there is ample experience with these techniques.^{2,4,8} Compared with histological biopsies, it is more easy to sample different tumor areas with fine-needle aspiration (FNA), and this may result in increased diagnostic accuracy, in particular when dealing with STBs with heterogeneous features on clinical imaging (CT and MRI).

Several different cell block methods can be used to process cell material thus collected.⁹ As an adjunct to routinely prepared smears or cell sediments, cytoblock techniques allow the application of IHC and molecular methods, expanding the diagnostic armamentarium. For this purpose, we and others^{8,10–14} have used the Cellient™ automated cell block system, by which cytotechnicians can make an automated cell

block within 1 hour, albeit with higher costs than that of traditional cell block techniques. As described previously, with the Cellient™ method, using methanol fixation instead of formalin, a broad array of diagnostically important antibodies can be applied to IHC after optimization of IHC protocols. In clinical cytology, the Cellient™ method has been used successfully for the characterization of tumor cells in serous fluids and FNAC material, for example, to characterize different carcinoma types or to diagnose metastatic melanoma.¹²

In this article, we report our first experience on the suitability of the Cellient™ method to diagnose several types of STB tumors, 8 primary lesions (5 of which were gastrointestinal stromal tumors) and 12 secondary recurrences or metastases, applying 9 diagnostically relevant antibodies that were not described in our earlier article of the Cellient™ method.

2 | MATERIAL AND METHODS

2.1 | Ethics statement

The study met the criteria of the code of conduct for responsible use of human tissue that is used in the Netherlands (Dutch federation of biomedical scientific societies; <http://www.federa.org>).

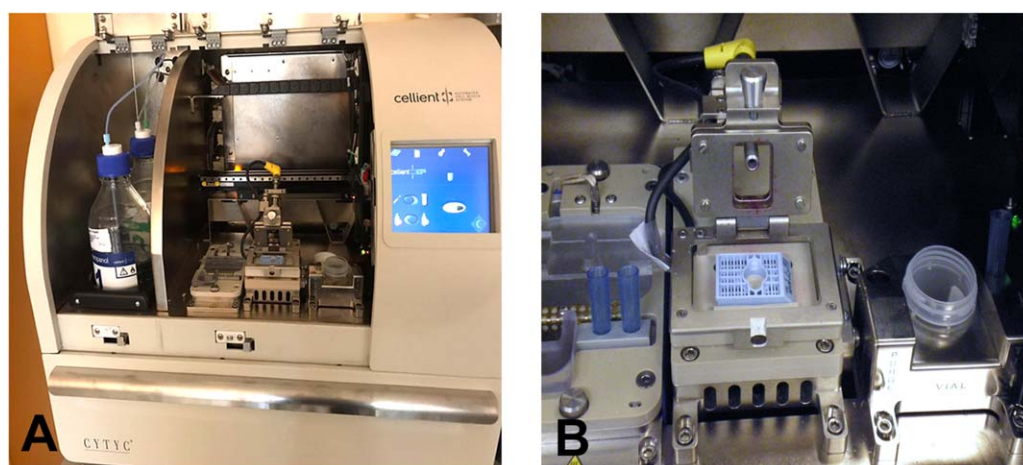


FIGURE 1 A closer look at the Cellient™ processor. Further details are found on the website: <http://www.hologic.com/products/clinical-diagnostics/instrument-systems/cellient-automated-cell-block-system> [Color figure can be viewed at wileyonlinelibrary.com]

2.2 | Cell samples

Cell samples of aspirations from soft tissue and bone tumors processed with the Cellient™ processor (Hologic, Marlborough, Massachusetts) between 2013 and 2016 were retrieved from the archives of the cytology laboratory of the pathology department of University Medical Center Groningen. Our cohort consisted of 20 consecutive cases, shown in Table 1, and included 12 EUS guided aspirations of deep-seated tumors (in the abdomen, retroperitoneal space, and mediastinum) and 9 FNAC specimens of superficial lesions. All but 1 EUS guided aspirations were performed with an EProCore needle (ECHO-HD-22-C; Echo Tip Ultra; Cook Medical, Bloomington, Indiana). In all EUS procedures and most FNAC aspirations, specimen cellularity had been checked on site by our cytotechnicians and in case of low cellularity, repeated aspirations were

done. The study cohort comprised 16 soft tissue tumor cases (5 were gastrointestinal stromal tumors) and 4 bone tumor cases (2 chondrosarcomas and 2 osteosarcomas). For evaluation of diagnostic performance, the cohort was divided in 2 groups, 8 primary lesions, and 12 secondary lesions (metastases or recurrences of tumors of which the histologic diagnosis was known). Histologic follow-up was available for all primary tumors, allowing correlation of cytological and histological diagnosis.

2.3 | Cellient™ cell block technique

Before being loaded into the Cellient™ processor (Figure. 1A,B), materials were washed in 1 mL Cytolyt™ Wash, centrifuged at 1000 g for 5 minutes, dissolved in 20 mL PreservCyt™ fluid and fixed for 20 minutes. One drop of the cell sediment was used to

TABLE 2 The 14 antibodies applied with Cellient™ cell block specimens

Antibody	Clone type	Clone	Manufacturer	Dilution	Pretreatment
Actin-SMA	Monoclonal	1A4	Ventana	R.T.U	No
Beta-catenin	Monoclonal	14	Ventana	R.T.U	CC1 52 min
Brachyury	Monoclonal	EPR18113	Abcam	1:400	CC1 36 min
Caldesmon	Monoclonal	h-CD	Dako	1:800	No
CD-117	Polyclonal	C-KIT	Dako	1:100	No
CD-34	Monoclonal	QBEND10	Ventana	R.T.U	CC1 92 min
CK-AE1/3	Monoclonal	AE1/AE3	Ventana	R.T.U	CC1 36 min + protease 4 min
Desmin	Monoclonal	DE-R-11	Ventana	R.T.U	CC1 64 min
DOG1	Monoclonal	SP 31	Ventana	R.T.U.	no
ER	Monoclonal	SP-1	Ventana	R.T.U.	No
HMB-45	Monoclonal	HMB45	Ventana	R.T.U	No
myf-4	Monoclonal	LO26	Monosan	1:25	CC1 64 min
S-100	Monoclonal	4C4.9	Ventana	R.T.U	No
SATB2	Monoclonal	4B10	Abcam	1:100	CC1 64 min

CC1, cell conditioning solution (Ventana), pre-treatment buffer, pH 8.4.

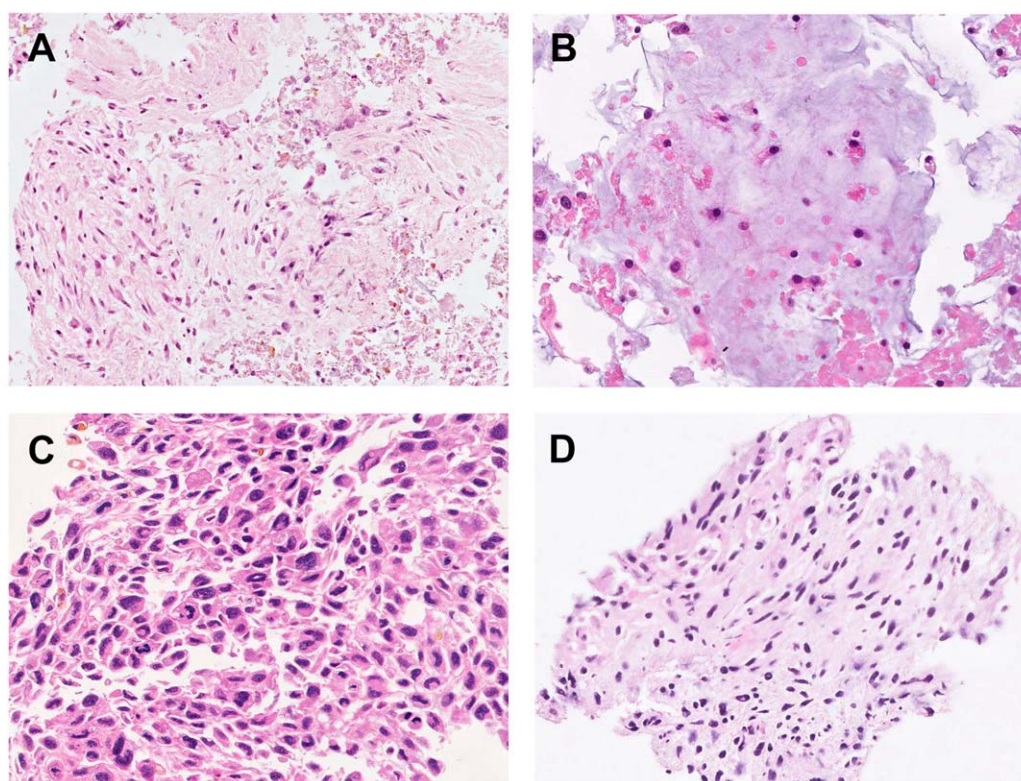


FIGURE 2 Microphotographs of H&E slides obtained from Cellient™ cell blocks. (A) collagen rich tissue with fibroblastic tumor cells in desmoid fibromatosis. (B) cartilaginous matrix with atypical hyperchromatic tumor cells in grade 2 chondrosarcoma. (C) pleomorphic tumor cells in pleomorphic undifferentiated sarcoma. (D) hyperchromatic spindled tumor cells in malignant peripheral nerve sheath tumor (original $\times 200$) [Color figure can be viewed at wileyonlinelibrary.com]

prepare a Giemsa-stained smear. In addition, 6 drops of the cell sediment were washed for 20 minutes in 1 mL Cytolyt Wash™ (Hologic), a low-dose methanol-based solution used to lyse erythrocytes and dissolve mucus. From this sediment, a Papanicolaou-stained microscopic thin layer slide was prepared with the ThinPrep T5000 processor. The remaining part of the cell suspension was rinsed twice in Cytolyt Wash™ solution and centrifuged again for 5 minutes at 1200 g, after which the pellet was fixed with Preserv-Cyt™ fluid for 20 minutes before the sample vial with Preserv-Cyt™ was put in the automated Cellient™ processor. The Cellient™ Automated Cell Block System is fully automated. It creates a paraffin-embedded cell block in <1 hour by means of a controlled vacuum that concentrates a layer of cells on a specially designed filter. Dehydrating and clearing reagents, including propanolol and xylene, are vacuum-drawn through the sample, which is subsequently embedded in paraffin and finished in an additional layer of paraffin; this makes it ready for histological sectioning. The vacuum-assisted filtration concentrates available cells within the final paraffin block. Eosin staining is used for visualization of the cell layer during sectioning. During sectioning of the Cellient™ cell blocks 10 paraffin sections of 4-mm thickness were prepared, and these were mounted on aminopropyltriethoxysilane (APES)-coated microscopic slides. One section was routinely stained with Hematoxylin and Eosin (H&E) for microscopic evaluation of specimen cellularity. The remaining unstained slides were available for IHC.

2.4 | Immunohistochemistry

The 14 antibodies (13 monoclonal, 1 polyclonal) applied in this study, including their commercial source, clone, and working dilution, as summarized in Table 2. Five diagnostically relevant antibodies had been evaluated in our earlier study of the Cellient™ method (CD117, AE1/3, ER, HMB-45, S-100). The 9 additional antibodies used to diagnose the STB tumors in this cohort were actin, beta-catenin, brachyury, caldesmon, CD-34, desmin, DOG-1, myf4, and SATB2. All IHC stains were performed in the Benchmark Ultra automated immunostainer (Ventana, Tuscon, Arizona) using the Ultraview detection system and validated by testing proper dilution of the antibody, need for CC1 antigen retrieval, and need for an 8 minutes amplification step in the IHC staining protocol, respectively (Table 2). All antibodies had been tested with at least 3 different Cellient™ cell blocks prepared from 3 different specimens. IHC results obtained with Cellient™ cell blocks were compared with IHC results obtained with corresponding formalin-fixed, paraffin-embedded (FFPE) tissue tumor material from the same patient as reference standard. Several antibodies required antigen retrieval with CC1 (cell conditioning buffer, pH 8.4) for optimal staining. CC1 with protease pretreatment proved to give the best results for cytokeratin antibody AE1–3. For all IHC staining, the Ventana Ultraview DAB detection kit was used with an amplification step of 8 minutes. Hematoxylin was used as a counterstain.

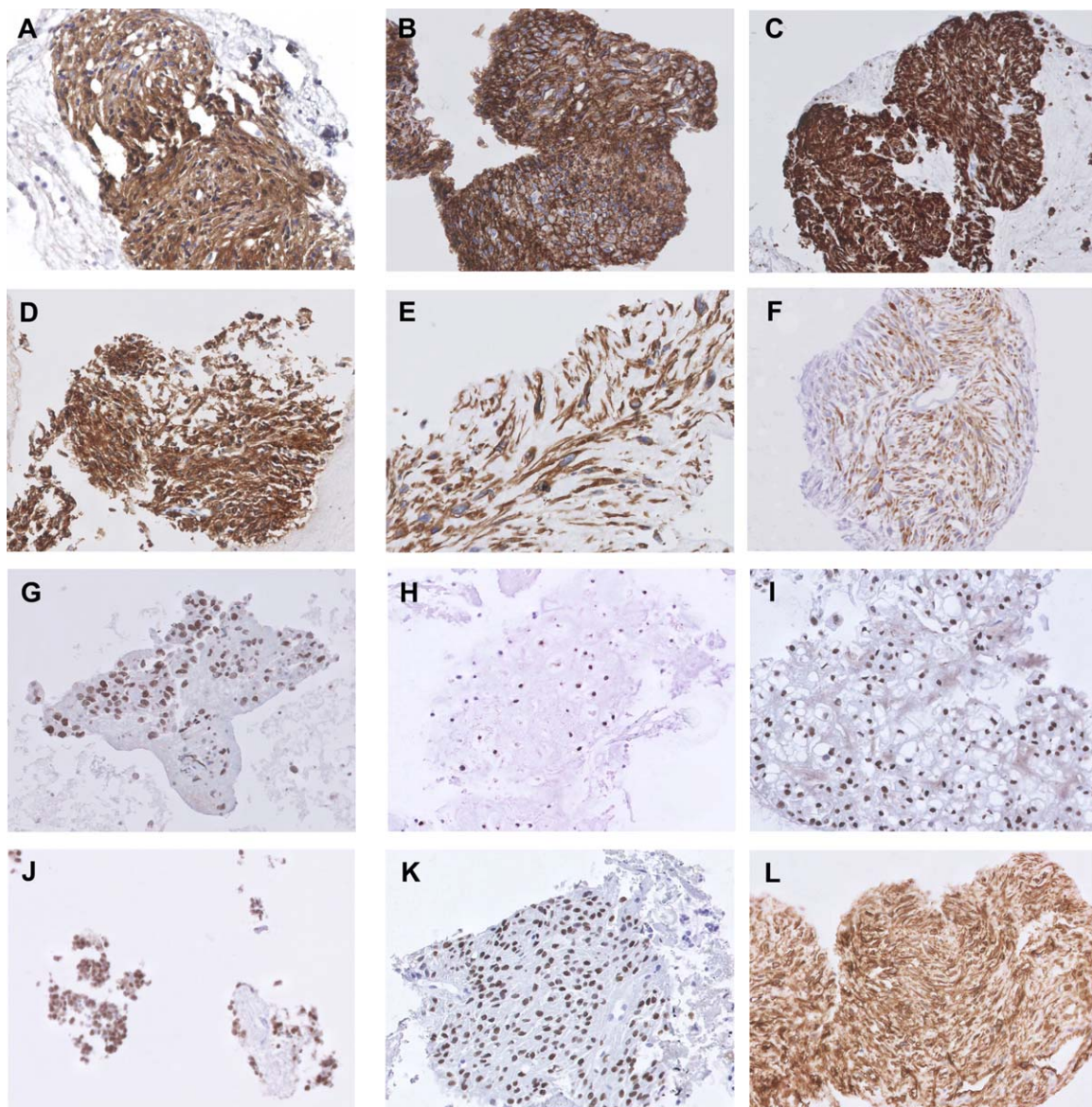


FIGURE 3 Microphotographs of IHC using 12 diagnostically relevant antibodies with the Cellient™ method. CD117 (A) and DOG-1 (B) in GIST. Desmin (C) and HMB-45 (D) in PEComa. SMA (E) and caldesmon (F) in leiomyosarcoma. Nuclear staining of osteosarcoma cells with SATB2 (G), S-100 (H) in chondrosarcoma, brachyury (I) in chordoma, myf4 (J) in alveolar rhabdomyosarcoma, and ER (K) in gynecologic leiomyosarcoma. CD34 (L) in dermatofibrosarcoma protuberans (original $\times 200$) [Color figure can be viewed at wileyonlinelibrary.com]

3 | RESULTS

In all 20 cases, the H&E stained sections of the Cellient™ material contained small tissue fragments. H&E histology of these small tissue fragments (microbiopsies) revealed that diagnostically relevant histological and cytological features could be examined properly, as shown in Figure 2. Fragments of desmoid fibromatosis consisted of collagen rich tissue with haphazardly arranged fibroblastic cells with round nuclei, nucleoli, and tapering eosinophilic cytoplasm (Figure 2A). Fragments of a grade 2 myxoid chondrosarcoma contained tumor cells with moderately atypical, hyperchromatic, single and double nuclei embedded in myxochondroid matrix (Figure 2B), whereas cellular fragments with pleomorphic and hyperchromatic tumor cells were encountered in recurrences of pleomorphic undifferentiated sarcoma and radiation

sarcoma (malignant peripheral nerve sheath tumor (MPNST)) (Figure 2C,D).

In the group of 8 primary tumors, we specifically diagnosed 5 spindle cell gastrointestinal stromal tumors (GISTs) by positive IHC for both CD117 and DOG-1. A PEComa of the rectum could be diagnosed after positive IHC for actin, desmin and HMB-45, and a mesenteric desmoid fibromatosis in a patient with familial adenomatous polyposis (FAP) syndrome showed focal nuclear positivity for beta-catenin, whereas a retroperitoneal leiomyosarcoma was diagnosed as it showed expression of the smooth muscle markers actin (SMA), desmin, and caldesmon. Thus, in all primary STB tumors, IHC performed on Cellient™ material provided clinically important information.

In the group of 12 secondary tumors of known STB primaries, we managed to confirm the presence of a local recurrence or metastasis in

TABLE 3 Comparison of commonly used cell block methods, as reviewed by Jain et al.⁹

Method	Advantage	Disadvantage	Utility	IHC	Molecular studies
Agar method	Inexpensive Better orientation of cell block	Inconvenient heat treatment process Heat-related artefacts possible, if not cooled as recommended	For any fluid or FNA	Optimum results for cytoplasmic and nuclear antigens	Suitable
Histogel method	Good cellular preservation and architecture	Tedious process as HistoGel needs to be converted and maintained in liquid state Possible heat-related artefacts	Useful in specimens with no visible sediment after centrifugation	Suitable	Suitable
Collodion bag method	Good cellular yield	Time-consuming preparation of bags Toxic ether fumes for storage	Friable tissues and fragments, specimens of scanty cellularity	Appropriate results	Suitable
Cellient method	Good cellular yield Uniformly distributed cells Improved cellular architecture and nuclear features Consistent results Automated method with reduced procedural time No cross contamination Minimal cell loss	Expensive machines and consumables Requires trained staff for cutting thin sections	Limited studies Useful in low-cellularity specimens Useful in cervical LBC	Good results with optimized IHC protocols, adjusted to methanol fixation (see Refs. 8 and 13)	High quality of DNA and RNA

Abbreviations: FNA, fine needle aspiration; IHC, immunohistochemistry; LBC, liquid-based cytology.

all cases. This group included 4 bone sarcomas. A metastasis of an osteosarcoma was positive for SATB2 and 2 secondary chondrosarcomas showed expression of S-100, whereas an inguinal lymph node metastasis of a sacral chordoma was confirmed by IHC for the transcription factor brachyury. The 8 secondary manifestations of soft tissue sarcomas included 2 radiation-induced sarcomas and 2 pleomorphic undifferentiated sarcomas (diagnosed solely on H&E morphology), 2 metastatic alveolar rhabdomyosarcomas (which were myf4 positive), a metastatic leiomyosarcoma of the uterus (in which actin and ER were positive), and a recurrent dermatofibrosarcoma protuberans (CD34 positive). Thus, in addition to H&E morphology, in 8 out of 12 cases, IHC on Cellient™ material provided incremental diagnostic information.

In both groups (primary and secondary tumors), IHC results in Cellient slides were concordant with those obtained in FFPE tumor biopsies or excisions/resections from the same patient.

Figure 3 depicts IHC results of all antibodies applied: CD117 (Figure 3A) and DOG-1 (Figure 3B) in GIST, desmin (Figure 3C) and HMB-45 (Figure 3D) in PEComa, SMA (Figure 3E) and caldesmon (Figure 3F) in leiomyosarcoma, SATB2 in osteosarcoma (Figure 3G), S-100 in chondrosarcoma (Figure 3H), brachyury in chordoma (Figure 3I), myf4 in alveolar rhabdomyosarcoma (Figure 3J), ER (Figure 3K) in metastatic gynecologic leiomyosarcoma, and CD34 (Figure 3L) in dermatofibrosarcoma protuberans.

4 | DISCUSSION

The Cellient™ is a fully automated device that produces a cell block within 1 hour based on a standardized protocol. This allows rapid diagnosis on the same day the specimen arrives in the lab instead of the following day, which is convenient in selected cases. A methanol-based

PreservCyt™ solution is used instead of formalin. Several research groups have mentioned that the cellularity of Cellient™ material is at least comparable to that in traditional cell blocks, whereas cytomorphological details, in particular chromatin structure, appear to be better.^{10,12,14} Advantages and disadvantages of commonly used cell block methods including Cellient™ have been amply reviewed by Jain et al.⁹ and are summarized in Table 3. We have noted that the Cellient™ cell blocks often contain small tissue fragments. In addition to H&E interpretation of histological features, IHC and molecular methods, for example, FISH or NGS, can be applied.⁸ In the cytology laboratory, immunostaining can be applied to cell smears, ThinPrep specimens, cytospin specimens, and cell blocks. In a UK NEQUAS quality control study, testing commonly used antibodies for a diagnosis of carcinoma, mesothelioma, melanoma, and lymphoma, it was found that the highest sensitivity was provided by cell blocks, followed by cytospin specimens, liquid-based cytology slides, and cell smears.¹⁵

Although cellular DNA and RNA are well preserved by methanol fixation, at the protein level, IHC protocols that are routinely used for FFPE material, have to be optimized and validated. We⁸ and Sauter et al.¹³ have extensively tested many different antibodies for Cellient™ material using the automated Ventana Benchmark immunostainer. In our initial study published in 2013, we showed that IHC performed on Cellient™ cell blocks could be applied to diagnostic algorithms that proved to be helpful in the discrimination of major tumor types (carcinoma, lymphoma, melanoma, and germ cell tumors), discrimination of carcinoma subtypes (adenocarcinoma, squamous-cell carcinoma, and neuroendocrine carcinoma), and determination of primary tumor site (eg, lung and breast) in cases of metastatic carcinoma. Notably, in a consecutive series of 100 cases, additional and clinically relevant information was obtained in 25% of serous fluid specimens and 29% of FNA specimens.⁸

To our knowledge, this is the first report on the use of Cellient™ cell blocks for the diagnosis of soft tissue and bone tumors. We stress that, in our and most sarcoma centers, cytology is only rarely applied for diagnosing tumors of soft tissue and bone. For the tumors in this series, EUS-FNA was tried to render a diagnosis of deep-seated primary tumors and FNA was used for superficially located metastatic or recurrent tumors of which the histologic diagnosis was known. We showed that by combining clinical presentation (including imaging studies), H&E morphology, and IHC, a diagnosis could be made in all 20 consecutive cases of tumors of soft tissue and bone. We evaluated 14 antibodies, 9 of which were not tested in our initial study. After optimization of factors influencing IHC results (in particular antigen retrieval conditions, amplification steps in the detection system kit) we managed to obtain excellent staining results for both cytoplasmic (eg, the smooth muscle markers, actin, desmin, and caldesmon) as well as nuclear antigens (eg, brachyury, myf4, and SATB2). In all 16 cases (all 8 primary tumors and 8 secondary tumors) where IHC was applied, a specific diagnosis could be made. For instance, brachyury, myf4, and SATB2, which are markers for notochordal, myogenic, and osteoblastic cell differentiation, respectively, allowed or confirmed a diagnosis of chordoma, alveolar rhabdomyosarcoma, and osteosarcoma.

Finally, cost considerations and budgetary constraints will determine the extent to which cytology laboratories use the rapid automated processing or more time-consuming traditional manual FFPE method to prepare cell blocks for H&E, IHC and/or FISH. Costs of the Cellient™ technique include purchase (50 000 US\$) and reagents (10 US\$ per specimen). Although the cost of the Cellient™ block technique is higher than that of a traditional cell block technique, we estimated that saved technician time is 30 minutes per specimen, using the time required to prepare an agar cell block as a reference standard. However, in our opinion, the cost of a new laboratory technique should be judged in the context of total cost of patient health care, including reduction of other diagnostic tests and patient life years saved, a cost analysis which is beyond the scope of this article.

In summary, we have shown that routine H&E staining and IHC of cell material processed with Cellient™ processor has the potential to accurately diagnose tumors of soft tissue and bone. In all 20 consecutive cases, important clinical information was provided, which translated into improved patient care. However, due to the small sample size, statistical analysis was not feasible, and a future study, testing appropriate antibodies on a larger number of cases, is needed to assess the real value of this method.

CONFLICT OF INTEREST

The authors have no conflict of interest.

DISCLOSURES

The authors made no disclosures.

ORCID

A. J. H. Suurmeijer  <http://orcid.org/0000-0003-1361-9454>

REFERENCES

- [1] Domanski HA. Fine-needle aspiration cytology of soft tissue lesions: diagnostic challenges. *Diagn Cytopathol*. 2007;35:768–773.
- [2] Kramer H, Sanders J, Post WJ, Groen HJ, Suurmeijer AJ. Analysis of cytological specimens from mediastinal lesions obtained by endoscopic ultrasound-guided fine-needle aspiration. *Cancer* 2006;108:206–211.
- [3] Layfield LJ, Baloch Z. Cytologic diagnosis of osseous lesions: a review with emphasis on the diagnosis of primary neoplasms of bone. *Diagn Cytopathol*. 2009;37:299–310.
- [4] van Hemel BM, Lamprou AA, Weersma R, Plukker JT, Suurmeijer AJ, van Dullemen HM. Procedure-related, false-positive cytology results during EUS-guided FNA in patients with esophageal cancer. *Gastrointest Endosc*. 2010;71:1130–1133.
- [5] Domanski HA, Akerman M, Carlén B, et al. Core-needle biopsy performed by the cytopathologist: a technique to complement fine-needle aspiration of soft tissue and bone lesions. *Cancer* 2005;105:229–239.
- [6] Kljanić J, Pierron G, Sastre-Garau X, Theocharis S. Value of combined cytology and molecular information in the diagnosis of soft tissue tumors. *Cancer Cytopathol*. 2015;123:141–151.
- [7] Layfield LJ, Armstrong K, Zaleski S, Eckardt J. Diagnostic accuracy and clinical utility of fine-needle aspiration cytology in the diagnosis of clinically primary bone lesions. *Diagn Cytopathol*. 1993;9:168–173.
- [8] van Hemel BM, Suurmeijer AJH, Aisner D. Effective application of the methanol-based PreservCyt™ fixative and the Cellient™ automated cell block processor to diagnostic cytopathology, immunocytochemistry, and molecular biology. *Diagn Cytopathol*. 2013;41:734–741.
- [9] Jain D, Mathur SR, Iyer VK. Cell blocks in cytopathology: a review of preparative methods, utility in diagnosis and role in ancillary studies. *Cytopathology* 2014;25:356–371.
- [10] Kruger AM, Stevens MW, Kerley KJ, Carter CD. Comparison of the Cellient™ automated cell block system and agar cell block method. *Cytopathology* 2014;25:381–388.
- [11] Montgomery E, Gao C, de Luca J, Bower J, Attwood K, Ylagan L. Validation of 31 of the most commonly used immunohistochemical antibodies in cytology prepared using the Cellient™ automated cell block system. *Diagn Cytopathol*. 2014;42:1024–1033.
- [12] Prendeville S, Brosnan T, Browne TJ, McCarthy J. Automated Cellient™ cytoblocks: better, stronger, faster? *Cytopathology* 2014;25:372–380.
- [13] Sauter JL, Grogg KL, Vrana JA, Law ME, Halvorson JL, Henry MR. Young investigator challenge: validation and optimization of immunohistochemistry protocols for use on Cellient cell block specimens. *Cancer Cytopathol*. 2016;124:89–100.
- [14] Wagner DG, Russell DK, Benson JM, Schneider AE, Hoda RS, Bonfiglio TA. Cellient automated cell block versus traditional cell block preparation: a comparison of morphologic features and immunohistochemical staining. *Diagn Cytopathol*. 2011;39:730–736.
- [15] Kirbis IS, Maxwell P, Fležar MS, Miller K, Ibrahim M. External quality control for immunocytochemistry on cytology samples: a review of UK NEQAS ICC (cytology module) results. *Cytopathology* 2011;22:230–237.

How to cite this article: Song W, van Hemel BM, Suurmeijer AJH. Suitability of the Cellient™ cell block method for diagnosing soft tissue and bone tumors. *Diagnostic Cytopathology*. 2018;46:299–305. <https://doi.org/10.1002/dc.23887>